

AFIAS Noro

- ✓ High accuracy
- ✓ Fast result
- ✓ Fully automated



AFIAS Noro

AFIAS Noro test is the high accuracy qualitative Noro virus test solution.

AFIAS Noro is a qualitative analytic reagent for the Noro viral antigen combination test based on high sensitivity fluorescence immunoassay (FIA).

Innovative technology enables AFIAS Noro accuracy of auto-analyzer, with easy and simple fully automated processing.

With AFIAS reader, AFIAS Noro test provides fast, simple automated test for Noro viral infection in stool sample

Accurate



Fast




Simple




AFIAS Noro

- ▶ **Improved sensitivity:**
High sensitivity fluorescence technology boasts 100 times of conventional fluorescence signal. It provides sensitivity of auto-analyzer. 4 times of ELISA and 16 times of rapid test.
- ▶ **Automated:**
Automated sample processing minimizes the operator exposure to sample.
- ▶ **Compact:**
Compact device can be used ubiquitously
- ▶ **Fast result:**
With 15 min of reaction time, test result right next to patient.

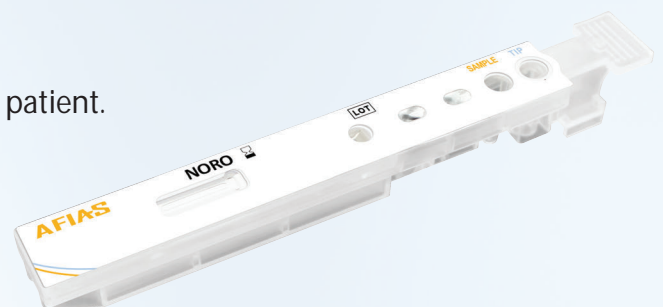

AFIAS Noro

		AFIAS Noro	ELISA	Rapid
Negative		-	-	-
Positive		+	+	+
Rbio Noro control	1	+	+	+
	1/10	+	+	+
	1/100	+	+	-
	1/200	+	+	-
	1/400	+	-	-
	1/800	+	-	-

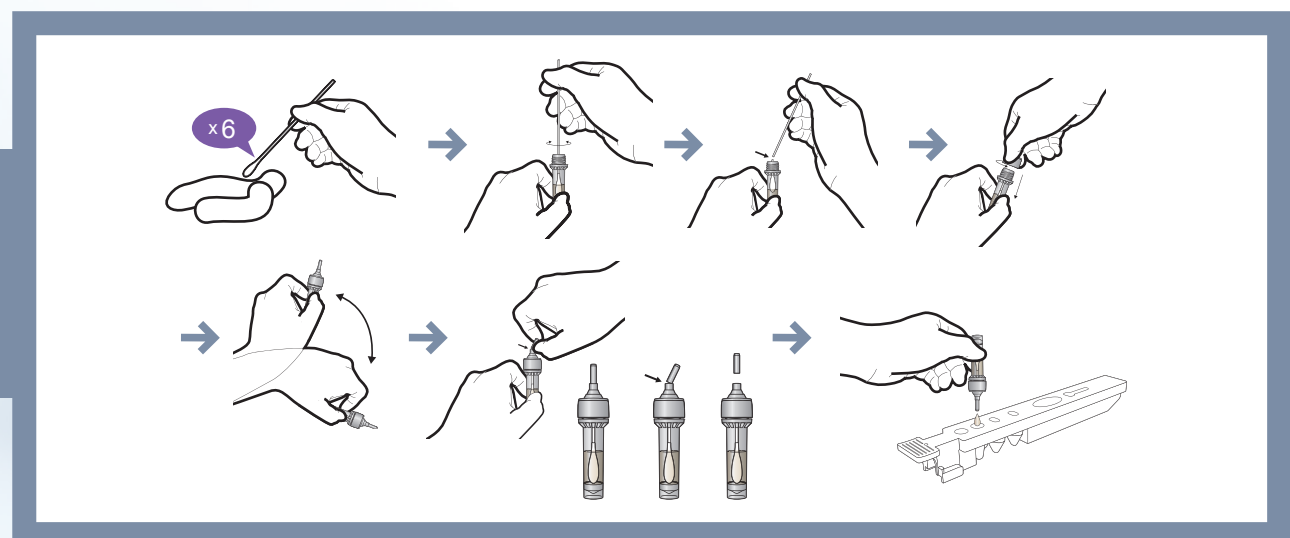

Item	Sensitivity	Specificity
AFIAS Noro	99%	87.6%

Item	AFIAS Noro
Cat #	SMFP-26
Test principle	FIA
Test type	Qualitative
Time	< 15 Min.
Sample	Stool
Storage	Cartridge: 2~8℃